



DIMAR

Medical Devices

DimAir Line

Ambulance
Emergency
& Hospital
Departments



EASY VENT MASK


EASY VENT MASK

“Easy Vent Mask” system for CPAP, recommended for the administration of air/oxygen mixtures in CPAP techniques at high flows on hypossiemic patients with Cardiogenic Aedema or Acute Respiratory Failure (ARF). It is provided with Venturi device, able to supply up to 120 Lt/min. and adjustable FiO₂ from 0,35 to 1.

Light, practical and quick, it can be employed in Ambulance, Emergency and hospital departments with only one oxygen source as distribution plant or cylinder.

TECHNICAL FEATURES

- Integrated Highly efficient single use Venturi
- Single patient
- Mask with anatomical inflatable cushion
- Seamless adjustable PEEP Valve from 5 to 20 cm/H₂O
- Possibility to administer FiO₂ from 0,35 to 1 (35% ÷ 100%)
- Integrated manometer, aimed to monitor and adjust the applied pressure continuously
- Connector for oxygen administration or EtCO₂ control
- 2-mt long tubes for oxygen supply
- Quick placement and fixation of the headgear
- Latex Free

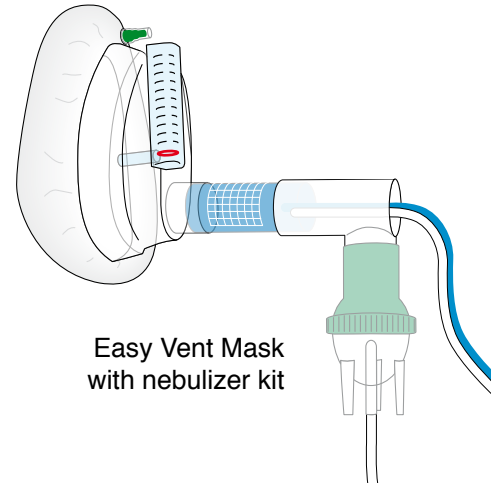

ADVANTAGES

- Extremely Quick, Practical and Efficient
- Continuous monitoring of the pressure applied even in case of leakages
- The efficiency of the system allows to administer pressures up to 20 cm/H₂O
- The reduced use of oxygen increases the cylinder autonomy
- Only one oxygen source is needed for its functioning, like distribution plant or cylinder
- Possibility of nebulization through specific kit (code 700/5065)

TABLE OF PERFORMANCES

Flow/Pressures/FiO₂
Obtained flow: about 60 Lt/min

Conc. O ₂ % approx.	40%	50%	60%
FLOWS O ₂ Lt/min.	A + B	A + B	A + B
PEEP cm/H ₂ O	5	7 10	6 18
	10	9 8	8 16
	15	11 6	10 14
OBTAINED FLOW approx.	60 Lt/min.		

REFERENCES

- Single pack
- Multiple pack with 10 pieces

CODE	SIZE
700/6412	PED/ADULT (S/M)
700/6411	ADULT (L)
700/6462	ADULT (XL)

Possibility to customize the different models, versions and packaging in response to customers' specific requirements.