



**DIMAR**

Medical Devices

100%  
MADE IN ITALY

# DimAir Line

**Ambulance  
Emergency  
Department**

5 POINTS  
HEADGEAR  
ATTACHMENT


LARGER  
NOSE ROOM

SINGLE PATIENT  
AND  
REPROCESSABLE  
VERSIONS

DEDICATED  
BRONCHOSCOPY  
PORT

DEDICATED  
NGT PORT

INTERCHANGEABLE ELBOWS  
VENTED/NON VENTED

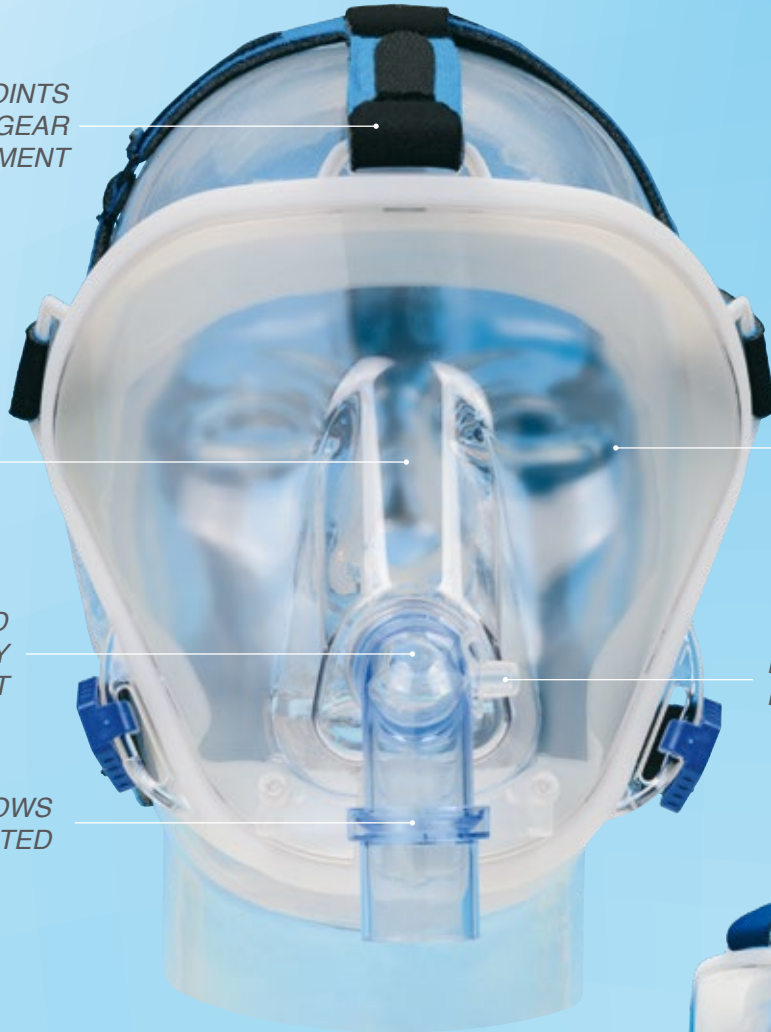

## TOTAL FACE DiMax

# DiMax TOTAL FACE MASK

DiMax is an innovative multi-purpose total face mask to perform Non Invasive Mechanical Ventilation in CPAP, BIPAP and PSV mode. Available in different sizes to fit paediatric and adult patient, DiMax can be used in combination with dual or single limb ventilators as well as CPAP high flow generators, simply changing the elbow connector. Unique features like 5 points headgear attachment, NGT airtight port, Bronchoscopy airtight port, extra nose room make DiMax the perfect interface for non invasive ventilation.

## 7 plus

- + attachments
- + comfort
- + sealing
- + practicality
- + versatility
- + safety
- + versions



The fifth point of attachment allows a better sealing onto the patient's forehead



Dedicated airtight NG tube port



Non-Vented bronchoscopy elbow

### TECHNICAL FEATURES

- 5 points of headgear attachment
- Larger nose room
- Dedicated ports for NG tube and Bronchoscopy probe
- Single patient and Reprocessable versions
- Interchangeable elbows:
  - 22F Non Vented (blue color)
  - 22M Vented (clear color)



Non-Vented "SE" elbow for dual limb ventilators



Vented "EE" elbow with safety valve for single limb ventilators

### ADVANTAGES

- Improved patient-ventilator synchronism
- Improved airtight performance
- Improved comfort regardless patient face contour or nose size
- Lesser headgear traction
- Lesser risk of skin lesions
- NG tubes and Bronchoscopy probes can be used without air leaks or interrupting the ventilation
- Reusable single patient version
- Reprocessable version sterilizable in autoclave @ 134° - Sterrad – ETO



Easy Vent Kit



Non-Vented elbow with safety valve

### REFERENCES

CODE*			SIZES
SINGLE PATIENT**	REPROCESSABLE***	REPROCESSABLE****	
700/6709	700/6739	700/6743	SMALL
700/6710	700/6740	700/6744	MEDIUM
700/6711	700/6741	700/6745	LARGE
700/6712	700/6742	700/6746	EXTRA LARGE

\* Including "Non Vented" blue elbow and "Vented" clear elbow

\*\* Single pack 5 pcs multi box

\*\*\* Single pack \*\*\*\* 5 headgears included