



Medical Devices

*DimAir* Line

# *The Helmet* **REVOLUTION**



Intensive Care  
Haematology  
Infectivology  
Home Care

## **ZIP HELMETS**



# ZIP HELMETS

Single patient Helmets for Non Invasive Mechanical Ventilation (NPPB), that can be used both in PSV-ASB, BIPAP and CPAP. Supplied with a pneumatic seal zip, they allow a quick and easy access to the patient. Comfortable and practical.


## TECHNICAL FEATURES

- 240° opening ZIP with pneumatic seal
- Made of biocompatible plastic material
- Additional fixing eyelets on the top of the helmet
- Symmetrical right and left connectors for an easy and practical connection with the equipment placed on any side
- Airtight access ports from 3 to 19 mm
- Straps made of soft and flexible material to reduce ulcerations
- Protective pads made of hygienic anti-bacterial CoolMax® material
- ISO 22M, 22Fconnectors
- Latex Free



## ADVANTAGES

- Quick access to the patient's face
- In PSV-ASB modality, the fixing « eyelets » improve the inspiratory trigger reaction and the synchronism between the patient and the ventilator
- Possibility to carry out cleaning, feeding and weaning operations without taking the helmet off
- The gas inflow and outflow can be inverted or connected from the same side ; in this way the patient can be placed on his side as well

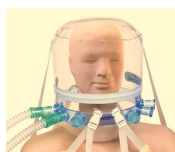

### NIMV 6R ZIP

Helmet for Non Invasive Mechanical Ventilation (NPPB), that can be used in PSV-ASB, BIPAP and CPAP.

6 ISO connectors placed on the rigid ring for a greater Efficiency, Practicality and Comfort. The connection to the rigid ring reduces the delay of the trigger reaction time <0,1"



### COMFORT ZIP

Helmet for CPAP. Devoid of the rigid ring, it prevents the formation of neck ulceration. The high comfort allows it to be used for prolonged periods both in hospital departments and In home care treatment.



### COMFORT NEO-PAED ZIP

CPAP Helmets for neonatal patients from 3 to 5 kg and paediatric patients from 5 to 15 kg. Devoid of the rigid ring, they prevent the formation of ulceration to the patient's neck

Fixing pad fabricated from medical approved anti-bacterial Coolmax® material.

NEONATAL	500/5447	3-5 Kg
PAEDIATRIC	500/5448	5-15 Kg

NIMV 6R ZIP	COMFORT CPAP ZIP	SIZE
600/5586	500/5678	XXS
600/5234	500/5441	XS
600/5235	500/5442	S
600/5236	500/5443	M
600/5237	500/5444	L
600/5238	500/5445	XL
600/5239	500/5446	XXL
600/5582	500/5568	3XL

- Single pack
- Multiple pack with 5 pieces

Possibility to realize different models, versions and packaging on customers' specific requirement.

2533GB/02



REG. N. MED 24039  
UNI EN ISO 9001:2008  
UNI EN ISO 13485:2012  
UNI EN ISO 14001 :2004



**DIMAR**

Via Galilei 6, 41036 Medolla -MO- I  
Tel: +39 0535 611336 Fax: +39 0535 611328  
Mail: dimar.srl@tiscali.it